


# Visible and Infrared Fourier Transform Spectroscopy (FTIR) • Beamline 1.4.2

Berkeley Lab • University of California

## Beamline Specifications

Photon Energy Range (cm <sup>-1</sup> )	Spectral Resolution (cm <sup>-1</sup> )	Spot Size (mm)	Availability
15–25,000 (0.002–3 eV)	0.125	1 (varies with coupling optics)	NOW

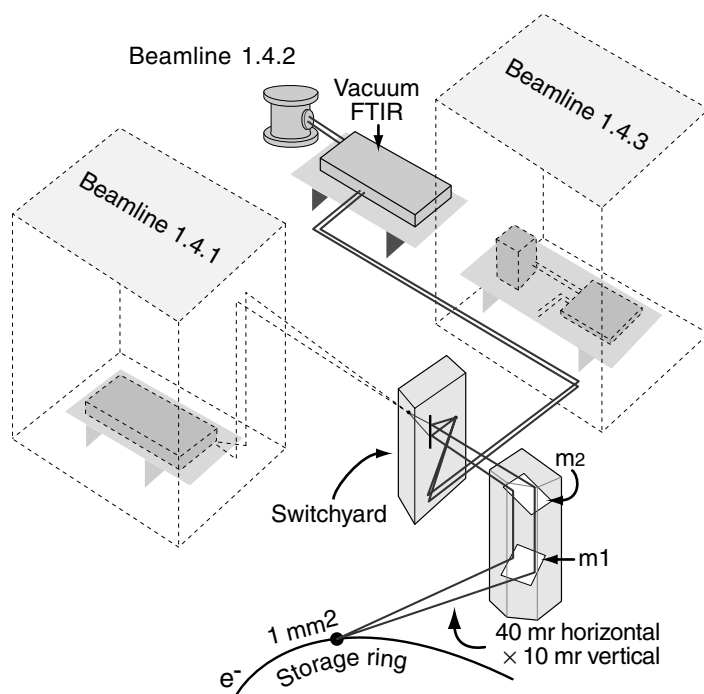

Schematic layout of Beamline 1.4.2.

**B**eamline 1.4 serves three distinct experimental stations designated as Beamlines 1.4.1, 1.4.2, and 1.4.3. Beamline 1.4.2 provides state-of-the-art Fourier transform infrared (FTIR) spectroscopy in the visible to far IR regions from 25,000 to 15 cm<sup>-1</sup>. All-reflecting optics focus the bend-magnet radiation into a switchyard. The switchyard houses optics that collimate the radiation and distribute it to the three experimental stations.


The synchrotron beam serves as an external source for a Bruker IFS66v/S vacuum FTIR bench with rapid-scan and step-scan capabilities. The beam, after being modulated by the moving mirror

in the Michaelson interferometer in the IFS66v/S, can be directed into an ultrahigh-vacuum (UHV) surface-science chamber. Any of several detectors placed on the opposite side of the surface-science chamber to measure the reflected signal may be operated remotely from the IFS66v/S. The IFS66v/S also has a standard internal sample area that may accommodate any standard FTIR sample accessory, including a LHe cryostat (2–475 K), attenuated total-internal-reflection cells (ATR), near-normal-incidence reflectivity adapters, solid pellets (KBr, polyethylene), and many other custom configurations.

The step-scan and fast-electronic capabilities of the IFS66v/S allow for time-resolved FTIR spectroscopy at time steps as short as 5 ns. In multi-bunch operation of the storage ring, pulses of light of approximately 44-ps duration illuminate the sample every 2 ns. Faster electronics than currently available but expected to be forthcoming should allow subnanosecond timing. A chopper is available for experiments at time scales longer than 10 ms.

Typical applications include pump-probe measurements (semiconductors, metastable states),

environmental science (adsorbates, bacteria, soil chemistry, remediation), biological materials (bioremediation, identification of biomolecules, tissue analysis), thin films and crystals (novel electronic materials and molecules), high-pressure systems (materials in diamond anvil cells), and measurements in high magnetic fields (reflectivity from high- $T_c$  materials). Visit [infrared.als.lbl.gov](http://infrared.als.lbl.gov) for more information. ■

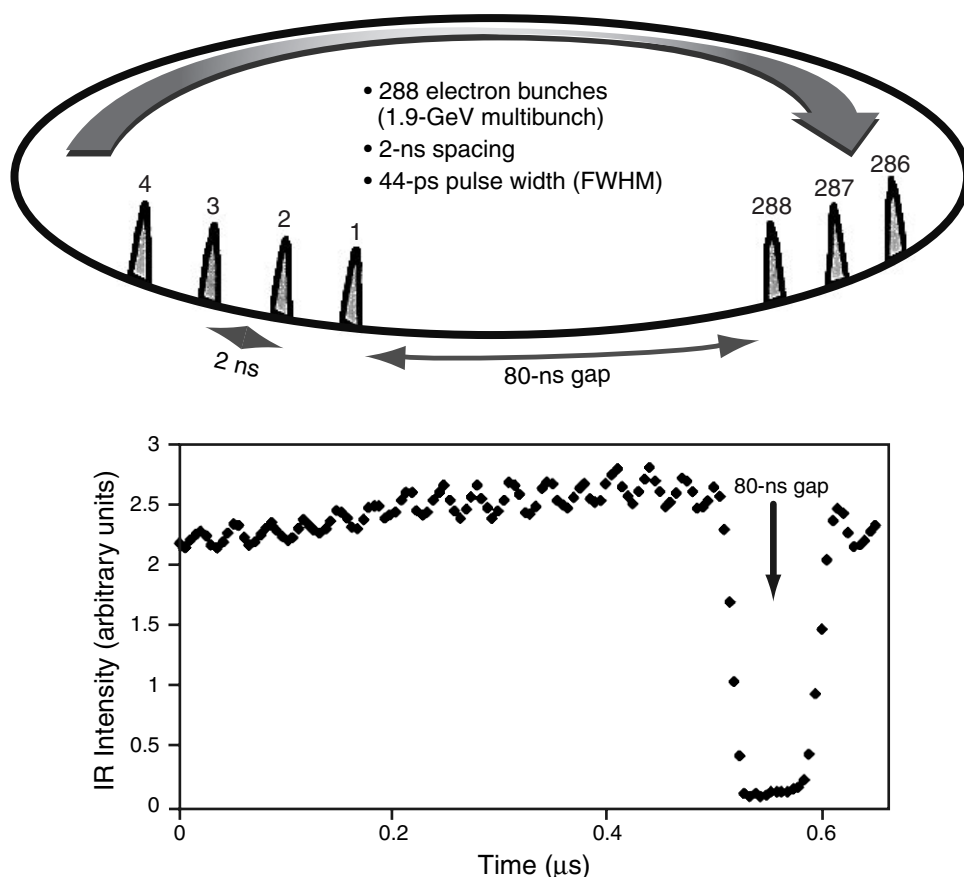

**Time-resolved IR spectroscopy on Beamline 1.4.2.** The step-scan and fast-electronic capabilities of the IFS66v/S allow for time-resolved FTIR spectroscopy at time steps as short as 5 ns. In the multibunch mode, the ALS produces pulses at intervals of 2 ns, except for a single long gap of 80 ns. This time-resolved plot of infrared intensity shows the 80-ns gap, but the time resolution of the current electronics does not fully separate the 2-ns spacing between pulses. Data courtesy of Mike Martin and Wayne McKinney (ALS).

**This beamline is available to independent investigators by submitting a proposal.**

#### For Beamline Information

Michael C. Martin  
Advanced Light Source  
Berkeley Lab, MS 7-222  
Berkeley, CA 94720  
Phone: (510) 495-2231  
Fax: (510) 486-7696  
Email: [mcmartin@lbl.gov](mailto:mcmartin@lbl.gov)

#### PRT Spokesperson

Wayne R. McKinney  
Advanced Light Source  
Berkeley Lab, MS 2-400  
Berkeley, CA 94720  
Phone: (510) 486-4395  
Fax: (510) 486-7696  
Email: [wrmckinney@lbl.gov](mailto:wrmckinney@lbl.gov)

#### To Obtain a Proposal Form

Bernie Dixon  
User Services Office Manager  
Advanced Light Source  
Berkeley Lab, MS 6-2100  
Berkeley, CA 94720  
Tel: (510) 486-6722  
Fax: (510) 486-4773  
Email: [alsuser@lbl.gov](mailto:alsuser@lbl.gov)