


High-Precision EUV Calibration and Standards • Beamline 6.3.2

Berkeley Lab • University of California

Beamline Specifications

Photon Energy Range (eV)	Photon Flux (photons/s/0.01%BW)	Spectral Resolution (E/ΔE)	Spot Size (μm)	Availability
50–1300	10 ¹¹ (at 100 eV)	7000	10 (v) × 300 (h)	NOW

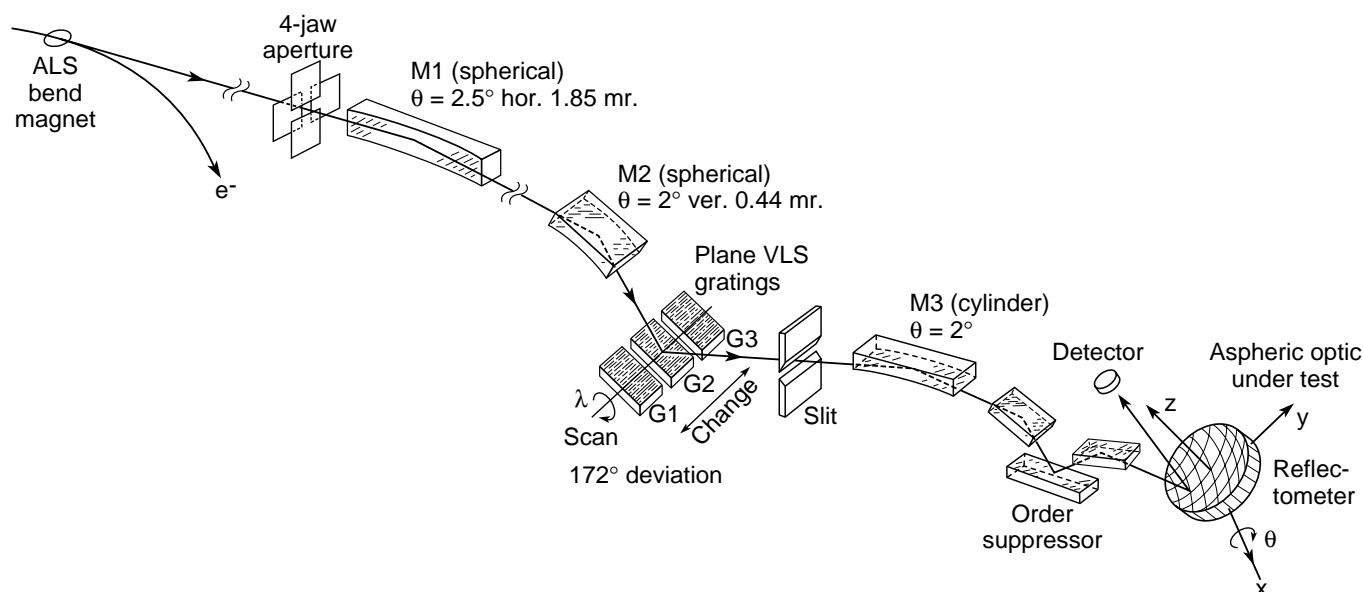

Schematic layout of Beamline 6.3.2.

Beamline 6.3.2 is a PRT-owned bend-magnet beamline dedicated to extreme ultraviolet (EUV) and soft x-ray reflectometry and scattering. The beamline is designed for high spectral purity and wavelength accuracy. Owned by the Berkeley Lab's Center for X-Ray Optics, the beamline is used for the characterization of optical components and reflective coatings for a variety of applications, including EUV lithography.


High spectral resolution is obtained using a variable-line-spaced plane-grating monochromator. The monochromator, designed and constructed by

the Center for X-Ray Optics, has a very compact design with three variable-line-spaced gratings, no entrance slit, and a fixed exit slit. Entrance-slitless operation is possible because of the small size and high stability of the ALS beam. A vertically deflecting spherical mirror converges the EUV beam onto the monochromator grating. Aberrations of the mirror are corrected by the varied line spacing of the grating, so that the spectral resolving power of 7000 is determined primarily by the small size of the ALS beam. The wavelength is scanned by simple rotation of the grating with respect to the

fixed monochromator exit slit. The light is focused onto the sample by the first horizontally deflecting mirror and a bendable refocusing mirror downstream from the monochromator. High spectral purity is achieved using a combination of filters and a triple-mirror "order-suppressor."

A permanent reflectometer experimental station is available. A two-circle goniometer permits the sample to be positioned in three dimensions to a precision of 1 micron and its angular position to be set to 0.002° . Samples of up to 20 cm in diameter

may be accommodated. An array of detectors including a photodiode, channeltron, and CCD camera are mounted on a rotating arm. Reflectometer performance specifications include high precision ($\Delta R/R \approx 0.2\%$) and high dynamic range (10^{10}). Space is available for additional experimental stations downstream from the reflectometer for a variety of applications, including spectroscopic studies in atomic, molecular, and materials science. Additional information is available at [www-cxro.lbl.gov/metrology](http://www.cxro.lbl.gov/metrology). ■

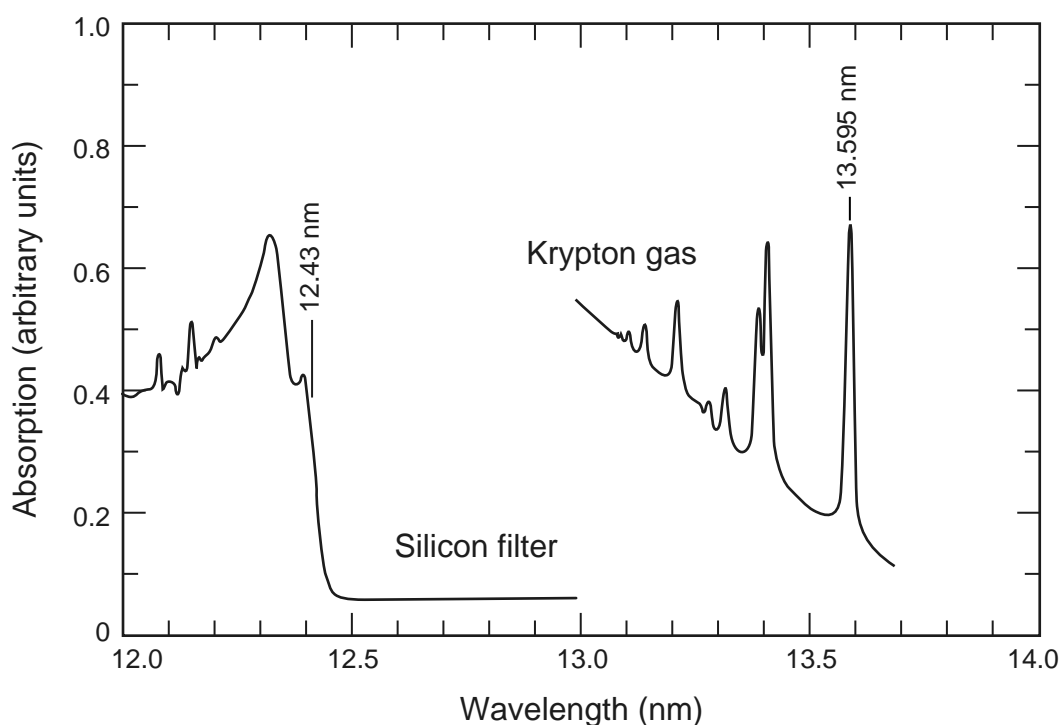

Absorption spectra of a silicon filter and of krypton gas used to calibrate the monochromator wavelength. Data courtesy of E. Gullikson (LBNL).

This beamline is available to independent investigators by submitting a proposal.

For Beamline Information

Eric Gullikson
Center for X-Ray Optics
Berkeley Lab, MS 2-400
Berkeley, CA 94720
Tel: (510) 486-6646
Fax: (510) 486-4550
Email: emgullikson@lbl.gov

Spokespersons

Stan Mrowka
Center for X-Ray Optics
Berkeley Lab, MS 2-400
Berkeley, CA 94720
Tel: (510) 486-4131
Fax: (510) 486-4550
Email: smrowka@lbl.gov

Ben Kaufmann
Center for X-Ray Optics
Berkeley Lab, MS 2-400
Berkeley, CA 94720
Tel: (510) 486-4079
Email: bbkaufmann@lbl.gov

To Obtain a Proposal Form

Bernie Dixon
User Services Office Manager
Advanced Light Source
Berkeley Lab, MS 6-2100
Berkeley, CA 94720
Tel: (510) 486-6722
Fax: (510) 486-4773
Email: alsuser@lbl.gov